


Anistel

SWIFT & FALCON

Surface sterilants for animal
healthcare environments



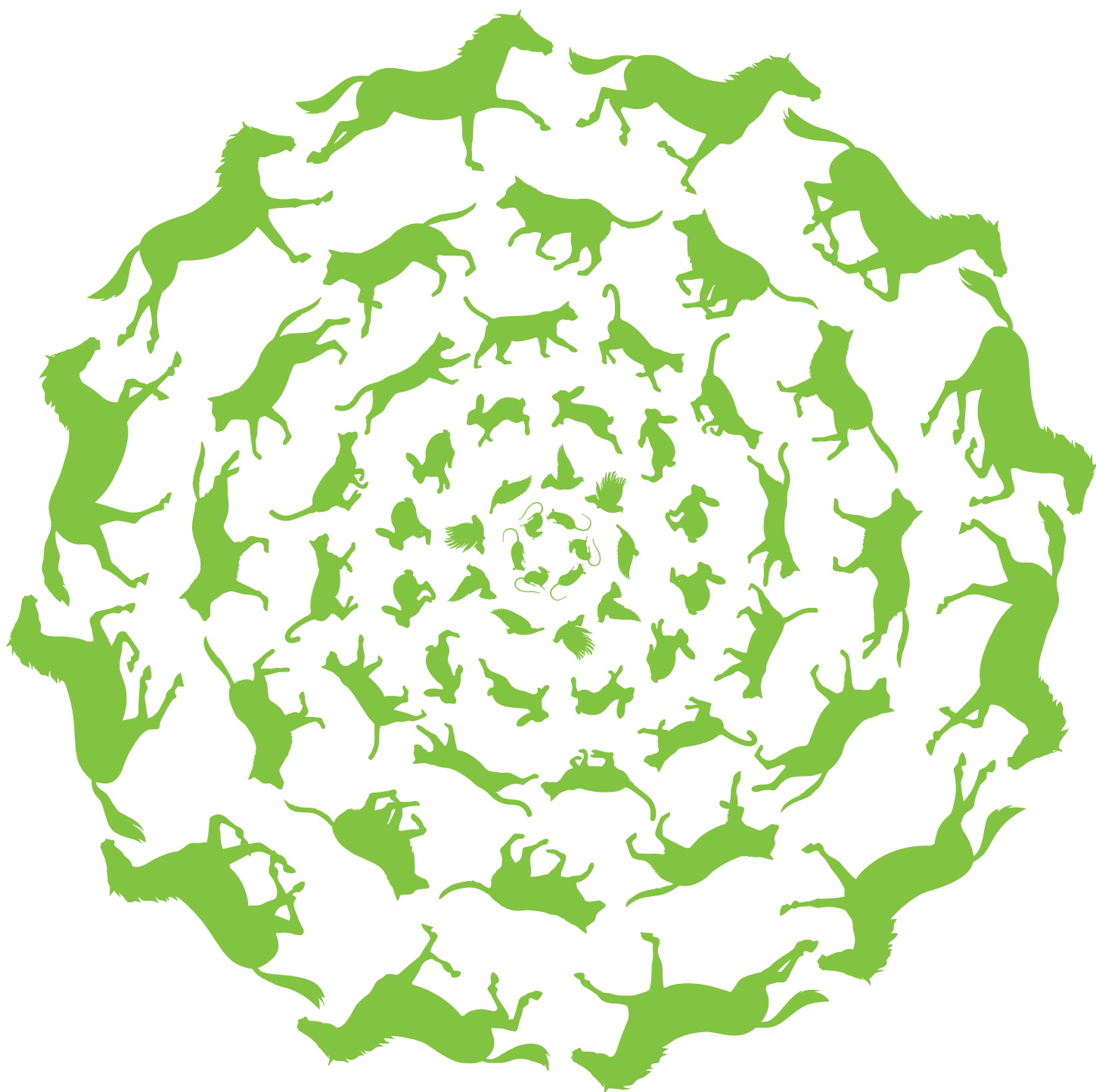

**Confidence for vets
and animals alike**

*Swift and Falcon utilise Tristel's
patented chlorine dioxide
chemistry, a well-documented,
highly effective and safe biocide*



Surface sterilants vets and animals can trust

Swift and Falcon are the latest additions to the Anistel range of Infection Control products. Both products are based on Tristel's patented chlorine dioxide (ClO₂) chemistry, which pioneered the replacement of glutaraldehyde in UK hospitals. Chlorine dioxide has been used by the NHS for 20 years and is a safe and powerful biocide that kills all microorganisms in very short contact times.

Why choose Swift and Falcon?

Classified and labelled in accordance with the Biocidal Products Directive 528/2012

Swift and Falcon are ideal for eliminating outbreaks and preventing closures:

- Effective against multi-drug resistant and spore forming organisms
- Practical and feasible contact times – 30 seconds for Falcon and 5 minutes for Swift
- Safe and easy-to-use

Suitable for the disinfection of hard surfaces

Swift and Falcon are designed for the rapid disinfection within:

- Operating theatres
- Isolation rooms
- Reception areas
- Consultation rooms
- Unoccupied animal habitats

Note: Swift and Falcon are not suitable for use on fabric or soft furnishings

Kill organisms in accordance with EN14885

Swift and Falcon are sporicidal, mycobactericidal, virucidal, fungicidal and bactericidal. They are effective against microorganisms of concern in animal healthcare, including:

- *Clostridium difficile* spores
- *Bacillus subtilis* spores
- MRSA
- Canine parvovirus
- Feline calicivirus (Norovirus)
- Foot-and-Mouth disease virus (FMDV)
- Swine Vesicular disease virus (SVDV)
- *Candida albicans*
- *Adenovirus*
- *Pseudomonas aeruginosa*
- *Staphylococcus aureus*
- *Escherichia coli*
- *Proteus vulgaris*
- *Aspergillus brasiliensis*
- *Mycobacterium avium*



SWIFT

Rapid action surface sterilant

Swift is a highly effective chlorine dioxide solution that is designed for the rapid disinfection of large hard surfaces within animal healthcare environments. When used as part of a daily or weekly deep clean regime, Swift is an excellent preventative method for outbreaks.

- Easy-to-use, safe and effective: one concentration kills all
- Sporicidal in 5 minutes
- For large hard surface areas
- No rinsing required

Chlorine dioxide's biocidal performance enables Swift to be sporicidal at one dilution in one short contact time of 5 minutes. This means that the user does not have to decide on different concentrations and dilution rates – one Swift sachet added to five litres of water will destroy all types of microorganisms.



Swift is available in boxes of 10 sachets

How to prepare and use Swift



Step 1

Pour 5 litres of cold water into container.

Step 2

Take one sachet.

Step 3

Fold sachet in half and squeeze one side to burst contents through centre seal. Contents will start to turn yellow.

Step 4

Swift is ready to use. Tear or cut sachet.

Step 5

Pour into five litres of water.

Step 6

Apply to surface with mop or cloth and leave to dry naturally. No rinsing is required.



FALCON

Highly effective surface sterilant spray

Falcon is a highly effective chlorine dioxide spray. It is designed for the rapid disinfection of small hard surfaces within animal healthcare environments. When used as part of a daily or weekly deep clean regime, Falcon is an excellent preventative method for outbreaks.




- Ready-to-use: no dilution required
- Sporicidal in 30 seconds
- For small hard surface areas
- No rinsing required

Falcon utilises a reusable trigger spray head for the distribution of chlorine dioxide. Falcon's base pack consists of one bottle containing sodium chlorite solution and one bottle containing a synergistic organic acid blend. This delivery format enables chlorine dioxide to be generated at point of use and at one concentration that kills all microorganisms.



Falcon is available in single packs

How to assemble and use Falcon

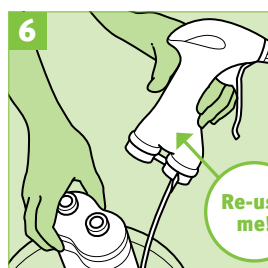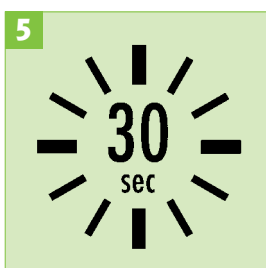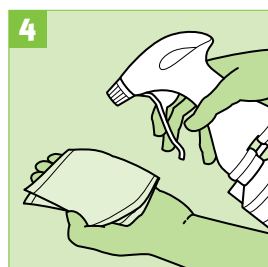

Step 1

Take one head and remove both bottle caps from the Falcon pack.

Step 2

Insert the the re-usable trigger head into the open bottle neck.

Step 3

Fasten the levers on each side. Falcon is now ready to be used.

Step 4

Apply Falcon to surfaces, avoiding ingress into equipment. Spread evenly with a wipe (Tristel Dry Wipes are recommended).

Step 5

Allow surface to dry naturally providing a contact time of at least 30 seconds. No rinsing is required.

Step 6

The head is designed to stop working when the refills are empty. Once the trigger head stops working, change the Falcon pack (even if there is liquid left in either side of the bottle). The trigger head can be re-used for the next Falcon pack.



Dry Wipes are available in tubs of 200



Follow us twitter.com/anisteluk
Join us facebook.com/anisteluk

Anistel is a good idea from **Tristel**

For more information call **+44 (0) 1638 721500** or email **mail@tristel.com**

www.tristel.com